PHILIPS

脉搏血氧仪
PULSE OXIMETER

操作指南 Quick Start Guide



型号(Model):F11TB/F16TB

感谢您选择飞利浦血氧仪，本页会快速帮助您使用本产品，如需了解更多请查阅说明书。
此产品商业代码为：
CO08F216/93、CO08F266/93

Thank you for choosing the Philips pulse oximeter. This page will quickly guide you through using this product. For more information,please refer to the Instructions for Use.
The commercial code for the product is:
CO08F216/93, CO08F266/93

扫码观看操作视频
Scan to watch the operation video

产品图例
PRODUCT LEGEND



APP连接
APP CONNECTION

01 下载并安装应用
DOWNLOAD AND INSTALL THE APP

扫描下方二维码或在应用市场搜索【飞利浦灵析】下载并安装APP
Scan the QR code below or search for [飞利浦灵析] in the APP store to download and install the APP

扫码下载APP
Scan the QR code to download the APP

如需不通过 APP 使用，详见背面说明
If you need to use it without the APP, please refer to the Instructions on the back


02 登录注册
REGISTRATION AND LOGIN

在APP主界面，根据向导完成账号登录，完善个人信息
On the main interface of the APP, complete account login according to the guide and complete personal information




03 添加设备
ADD DEVICE

打开蓝牙，在APP内点击【我的】→【添加设备】进入搜索设备界面
Turn on Bluetooth, tap [我的] → [添加设备] in the APP to enter the device search screen

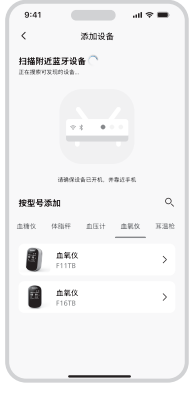

04 准备设备
PREPARE DEVICE

捏开血氧仪或按下按键，血氧仪开机
Pinch open the oximeter or press the button to power on

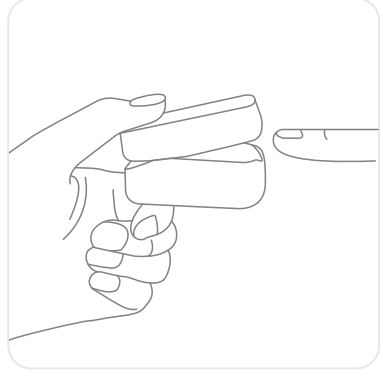

05 连接设备
CONNECT DEVICE

查找到需要绑定的设备进行连接
Search for the device to pair and connect



06 APP数据回溯
APP DATA BACKTRACKING


进入APP查看数据分析和历史数据
Enter the APP to view data analysis and historical data



仪器操作指南
INSTRUMENT OPERATION GUIDE


01 设备开机
POWER ON

捏开血氧仪或按下按键，血氧仪开机
Pinch open the oximeter or press the button to power on

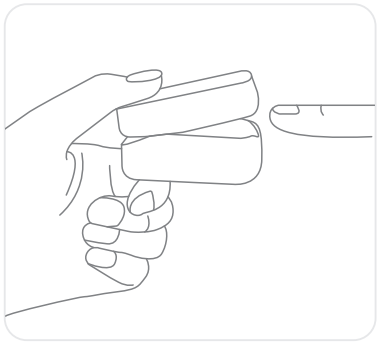


02 佩戴设备
FIT DEVICE

伸入手指前，请确保手部清洁，手指需完全伸入血氧仪腔体后（长指甲或涂有指甲油等会影响检测），即可使用设备
Before inserting your finger, please ensure your hands are clean. Fully insert the finger into the oximeter cavity before using the device. Long nails or nail polish may affect the measurement.

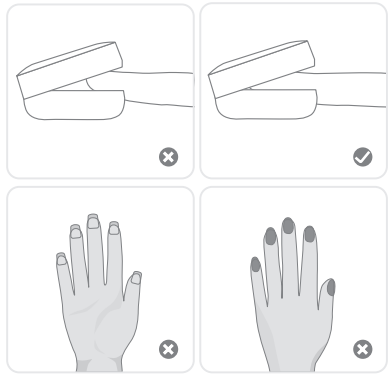


03 等待测量
WAIT FOR MEASUREMENT

等待3~6s机器识别，识别时血氧和脉率呈直线
Wait 3~6 seconds for the machine to recognize. During recognition, the blood oxygen and pulse rate will display as a straight line

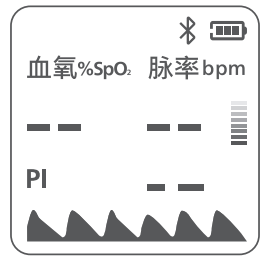

04 数据测量
DATA MEASUREMENT

屏幕开始呈现数字即开始测量，此过程为持续性过程，手指拔出即结束本次测量
Measurement starts when numbers appear on the screen. This is a continuous process; removing the finger ends the current measurement

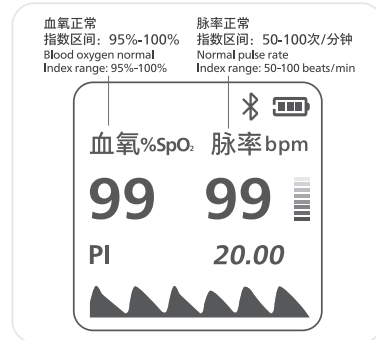

备注：此标准为临床医学与健康监测领域通用的标准参考范围
Note: This standard is a general reference range used in clinical medicine and health monitoring

05 数据回溯
VIEW HISTORICAL DATA

双击按键进入历史数据页面，单击按键进行跳行和翻页
Double press the button to enter the history data page. Press the button to move between rows and pages

历史数据		
	血氧(平均)	脉率(平均)
数据1:	97	91
数据2:	97	91
数据3:	35	250
数据4:	97	91
退出		


06 设备设置
SETTING DEVICE

长按按键进入设置页面，单击按键进行跳行和翻页，长按按键进行选项切换，最后长按【退出】键进行参数保存和退出
Press and hold the button to enter the settings page. Press the button to move between rows and pages, press and hold to switch options, and finally press and hold the [Exit] to save parameters and exit

Setting 设置	
Brightness 亮度	2
Alarm 报警	off 关闭
Monitoring 监测	off 关闭
Language 语言	en 英文
TH Setting	阈值设置
Restore	恢复默认设置
EXIT 退出	

07 挂绳安装
LANYARD INSTALLATION

根据图例将挂绳安装到设备上
Attach the lanyard to the device as shown in the device

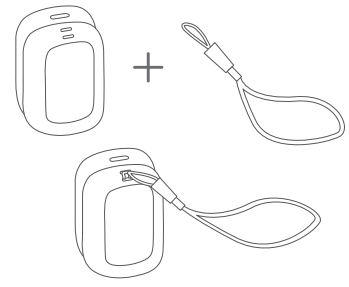

08 注意事项
PRECAUTIONS

请使用Type-C充电线搭配5V⎓1A(5W)的适配器进行充电，本产品不适用于大于20W的快充充电头
Please use a Type-C charging cable with a 5V⎓1A (5W) adapter for charging. This product is not suitable for fast chargers greater than 20W

